ALLMOBILITY TRADING S.r.l.

"BUFFALO 150" BOARD

USER TRANSFER SUPPORT



CE

Product complying with the requirements of Directive 2007/47/EC.

INSTRUCTIONS FOR USE
Version 001-16

"Kid-Man" is legitimate trademark of "Kasko Group" Ltd. All rights reserved. ©2025 "Kasko Group" Ltd.

INDEX

1. PRODUCT DESCRIPTION 3
2. TECHNICAL DATA 3
3. GENERAL INFORMATION 4
4. INSTRUCTIONS FOR USE 5
5. DISPOSAL 7
6. WARRANTY - ASSISTANCE 7

"Kid-Man" is legitimate trademark of "Kasko Group" Ltd. All rights reserved. ©2025 "Kasko Group" Ltd.

1. PRODUCT DESCRIPTION

The product described in these instructions is an easy-glide synthetic fibre board for user transfer.

It is designed to be easily transportable, in fact it is only 4 mm thick and has a handle hole.


The *BUFFALO 150* board is brightly coloured and has a double notch for greater stability and to ensure maximum safety during transfer.

The composite fibre ensures high resistance and avoids harmful chipping of the surface which occurs over time with other materials.

The *BUFFALO 150* board is an aid characterised by flexibility and versatility of use.

It can be easily used for small transfers in any indoor environment, for example in the bedroom, living room, kitchen or bathroom; or outdoors, for example in a car, on a bench or chair, for moving from or to the wheelchair.

Some illustrations:

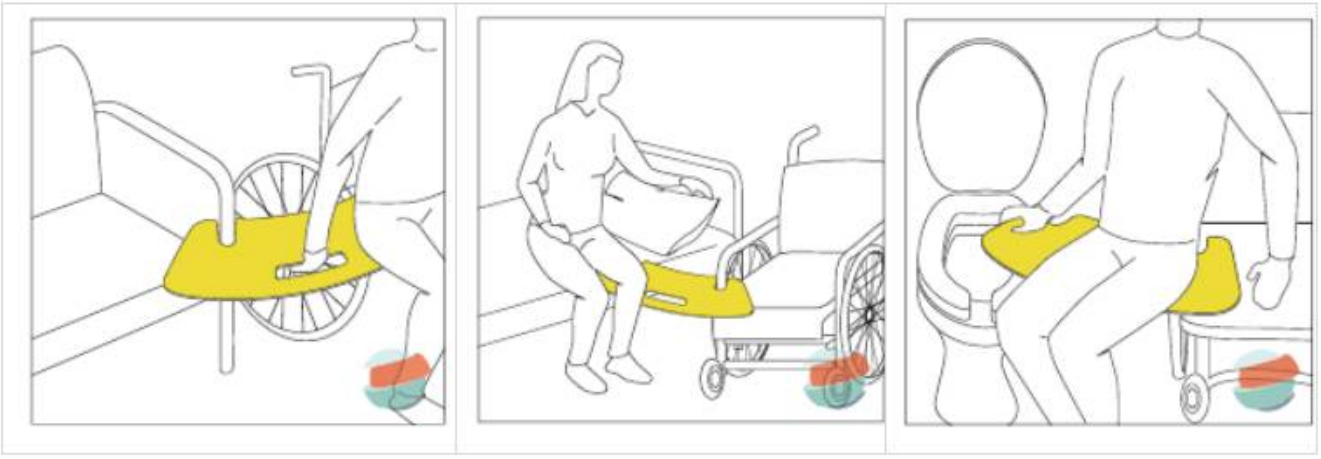



2. TECHNICAL DATA

"Kid-Man" is legitimate trademark of "Kasko Group" Ltd. All rights reserved. ©2025 "Kasko Group" Ltd.

"BUFFALO 150" BOARD

DIMENSIONS	WEIGHT (Kg)	BEARING CAPACITY (Kg)
 *LENGTH = 77cm* *WIDTH = 32.5cm* *THICKNESS = 0.4cm*	*1.23*	*150* *distributed over a surface of 40 cm, with resting weight of 25 kg on the sides*

STORAGE TEMPERATURE
ADMISSIBLE TEMPERATURE RANGE: *MINIMUM = -5 °C* *MAXIMUM = +40 °C*

3. GENERAL INFORMATION

"Kid-Man" is legitimate trademark of "Kasko Group" Ltd. All rights reserved. ©2025 "Kasko Group" Ltd.

The board has been designed and made to be safe in the use for which it is intended, as long as it is used and stored as indicated in these instructions.

Given below are the main instructions and correct storage by the user to as to reduce situations of risk.

 = PAY ATTENTION!

 = FORBIDDEN ACTION

	**The board is a medical aid for facilitating the transfer of a disabled user from a wheelchair to another seat.** **DO NOT USE THE BOARD FOR ANY PURPOSE OTHER THAN THE ONE DESCRIBED IN CHAPTER 2.**
	**DO NOT USE THE BOARD IF IT SHOWS ANY SIGN OF DAMAGE, DEFORMATION, CRACKS, THINNING, BREAKAGES.**
	**DO NOT USE THE BOARD TO TRANSFER WEIGHTS/PERSONS HEAVIER THAN THE MAXIMUM ALLOWED BEARING CAPACITY (150 Kg).**
	**DO NOT USE CORROSIVE PRODUCTS TO CLEAN THE BOARD.** **ALL IT NEEDS IS A DAMP CLOTH, SURGICAL SPIRIT OR ORDINARY SOAP.**
	**FOLLOW THE INSTRUCTIONS FOR USE PROVIDED ON THE BOARD AND IN CHAPTER 4.**
	**MAKE SURE THAT THE DISTANCE TO THE POINT OF ARRIVAL IS IDEAL FOR A CORRECT AND STEADY POSITIONING OF THE BOARD.**
	**WHEN POSSIBLE, AS WELL AS ONTO THE ARM OF THE WHEELCHAIR, HOOK THE BOARD ALSO ONTO A SUPPORT ON THE POINT OF ARRIVAL.**
	**IF THE INSTRUCTIONS ON THE LABEL APPLIED ON THE BOARD ARE SO WORN THAT THEY ARE NO LONGER LEGIBLE, ASK THE MANUFACTURER FOR A NEW COPY.**

4. INSTRUCTIONS FOR USE

"Kid-Man" is legitimate trademark of "Kasko Group" Ltd. All rights reserved. ©2025 "Kasko Group" Ltd.

The *BUFFALO 150* board is a support for transferring the user from the wheelchair to another position, for example to the bed, couch, car, etc.

It does not have any packaging, so it is bought, used, kept and carried just as it is.

Shown below are some of the most frequent uses and how to perform them, as shown also on the label applied on the accessory.



TRANSFER TO THE CAR

Get as close as possible to the seat, secure the board to the arm of the wheelchair and ensure that it is resting firmly on the seat.



TRANSFER TO THE BED

Get as close as possible to the bed, secure the board to the arm of the wheelchair and ensure that it is resting firmly on the mattress.



TRANSFER TO A CHAIR/COUCH

Get as close as possible to the chair/couch, secure the board to the arm of the wheelchair and to the armrest of the chair/couch (if possible) and ensure that it is resting firmly on the seat.



TRANSFER TO THE TOILET

Get as close as possible to the toilet, secure the board to the arm of the wheelchair and ensure that it is resting firmly on the seat.

ALWAYS FOLLOW THE INSTRUCTIONS IN CHAPTER 3.

"Kid-Man" is legitimate trademark of "Kasko Group" Ltd. All rights reserved. ©2025 "Kasko Group" Ltd.

5. DISPOSAL

When the board is no longer to be used, dispose of it correctly without dumping it in the environment.

Follow the rules for disposal in force in the place of use of the accessory (ask the local authorised disposal companies).

6. WARRANTY - ASSISTANCE

Listed below are the general terms of the warranty applied by the Manufacturer **ALLMOBILITY TRADING S.r.l.** .

The warranty applies to all manufacturing and material defects, duly confirmed at the Company's discretion.

***The Warranty does not cover** damage caused by:*

- *Normal wear and tear*
- *Deterioration of parts with purely aesthetic functions*
- *Repairs or modifications carried out by unauthorised persons or centres*
- *Use of the product not complying with the warnings given in this manual*
- *Damage caused by unsuitability of the environment in which the product is used and by phenomena not depending on its normal operation.*

In case of sale the warranty has a duration of 1 (one) year.

For any problem or to request information, contact the Manufacturer at the addresses on the first page.

"Kid-Man" is legitimate trademark of "Kasko Group" Ltd. All rights reserved. ©2025 "Kasko Group" Ltd.

7. VERSIONE BUFFALO GRIP 150


The BUFFALO GRIP 150 version (code All03 grip) **features two high-grip areas on the bottom.**

This stabilises the end surfaces of the board when they are in contact with soft parts, such as the padding or fabric of a chair, armchair, car or wheelchair cushion.

This allows some types of users to make the transfer according to their own preference.

The **material** is a "cushion grip" H3408 s, made of 51% polyurethane and 49% polyester, with an adhesive acrylic base, 1 mm thick.

The **shape and dimensions** (in mm) are those shown in the figure.

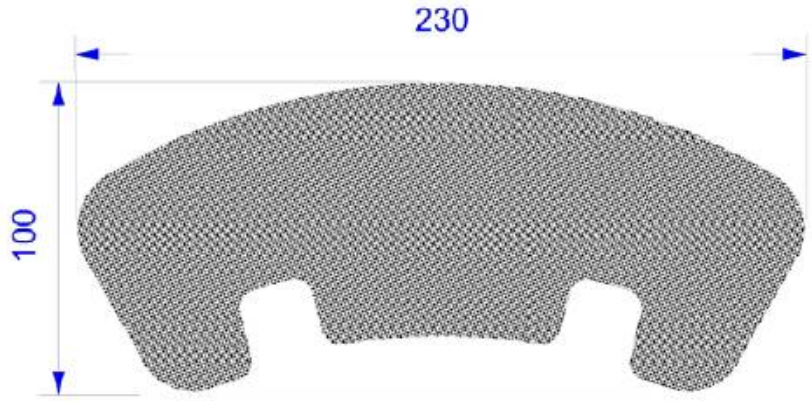

Do not use the board with the grip at the top (in contact with the body) during transfers.

"Kid-Man" is legitimate trademark of "Kasko Group" Ltd. All rights reserved. ©2025 "Kasko Group" Ltd.